



Edison Medical

UniBase™ System

By Nikola Khoury

Multi-Clicq™ System

ZIMAX® WP

Internal Hex



Digital Dentistry

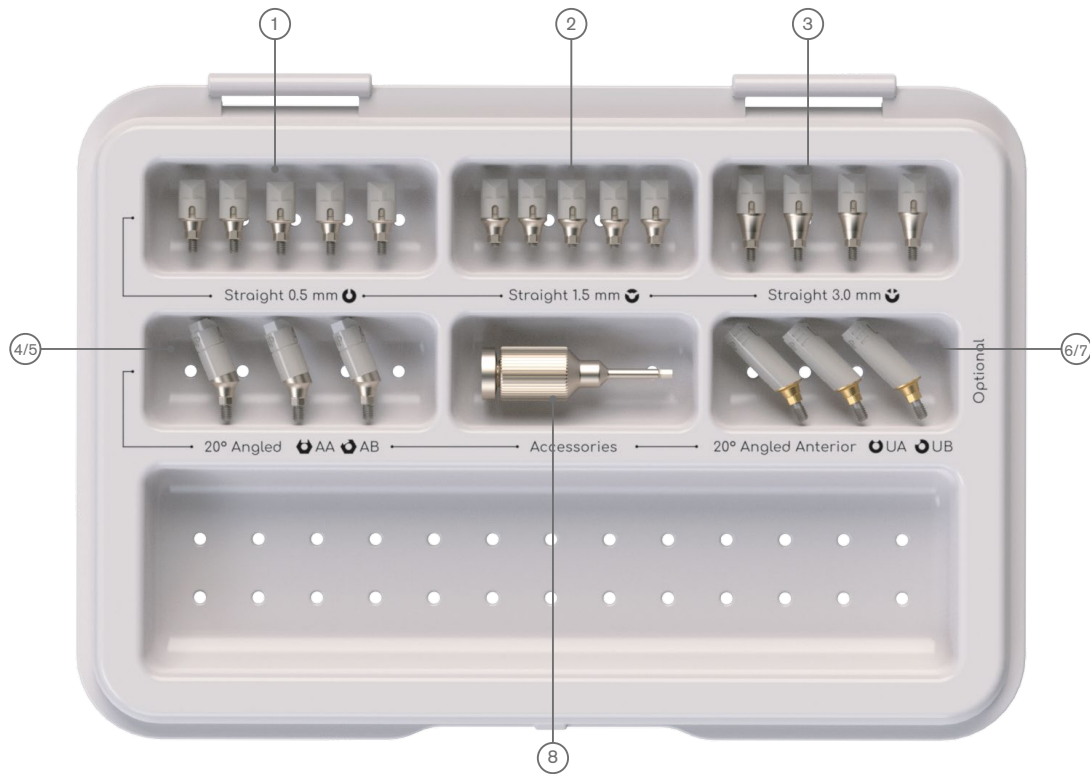
zinon

UniBase™ System

By Nikola Khoury

Abutment Level Coded Scan Posts

Patent Pending



Scan Post - Professional & Starter Kit

Internal Hex	Basic Kit	Starter Kit	Professional Kit
Zimax® WP	SP-SV-WPBK	SP-SV-WPSK	SP-SV-WPPK

Scan Post - Basic Kit

- 1 - Straight Scan Post H 0.5mm - 5 pcs
- 2 - Straight Scan Post H 1.5mm - 5 pcs
- 3 - Straight Scan Post H 3.0mm - 4 pcs
- 4 - Hand Driver



Scan Post - Professional Kit

- 1 - Straight Scan Post H 0.5mm - 5 pcs
- 2 - Straight Scan Post H 1.5mm - 5 pcs
- 3 - Straight Scan Post H 3.0mm - 4 pcs
- 4 - 20° Angled Scan Post - Standard AA Type - 2 pcs
- 5 - 20° Angled Scan Post - Standard AB Type - 1 pcs
- 6 - 20° Angled Scan Post - Anterior UA Type - 2 pcs
- 7 - 20° Angled Scan Post - Anterior UB Type - 1 pcs
- 8 - Hand Driver

Scan Post - Starter Kit

- 1 - Straight Scan Post H 0.5mm - 3 pcs
- 2 - Straight Scan Post H 1.5mm - 3 pcs
- 3 - Straight Scan Post H 3.0mm - 2 pcs
- 4 - 20° Angled Scan Post - Standard AA Type - 1 pcs
- 5 - 20° Angled Scan Post - Standard AB Type - 1 pcs
- 6 - 20° Angled Scan Post - Anterior UA Type - 1 pcs
- 7 - 20° Angled Scan Post - Anterior UB Type - 1 pcs
- 8 - Hand Driver

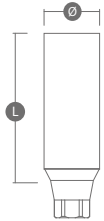
UniBase™ System

By Nikola Khoury

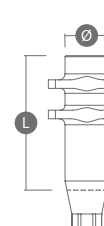
Digital and Manual Impression

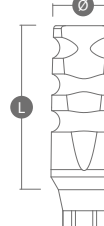


Digital Scan Body Classic Implant Level Scan

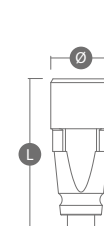
	
L	9.3mm
Ø	5.2mm
Internal Hex	
Zimax® WP	SB-SV-WP

Open Tray Transfer

	
L	12.5mm
Ø	4.5mm
Internal Hex	
Zimax® WP	OT-SV-WPN

		
L	9mm	12mm
Ø	3.8mm	4.5mm
Internal Hex		
Zimax® WP	CT-SV-WPSN	CT-SV-WPLN

Digital Implant Analog

	
L	10mm
Ø	4.5mm
Internal Hex	
Zimax® WP	DIA-SV-WP

UniBase™ System

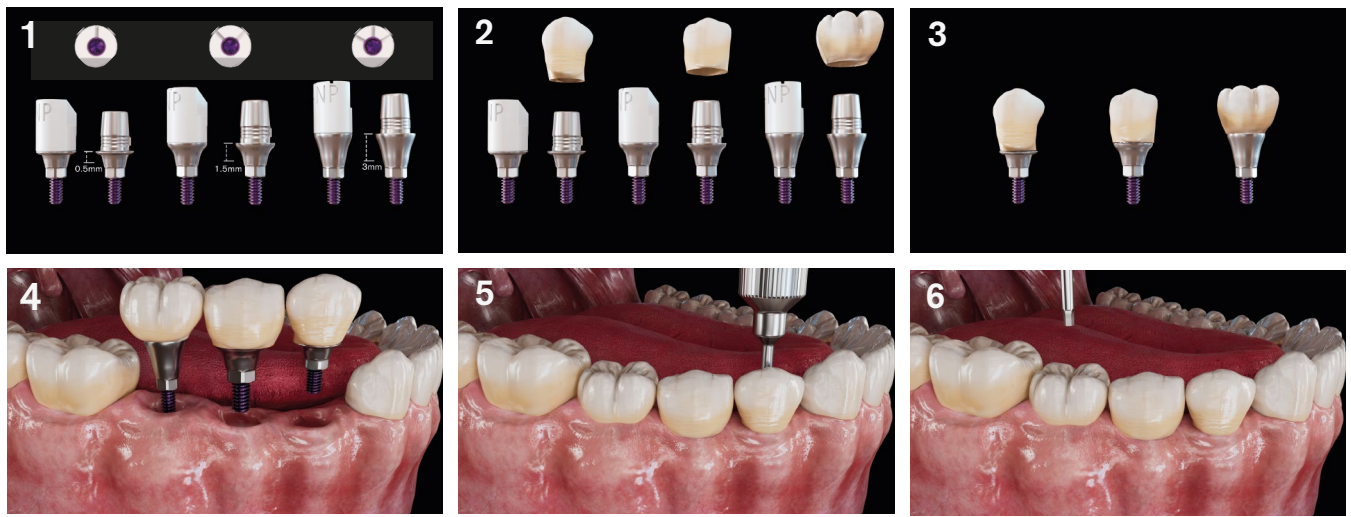
By Nikola Khoury

Straight Titanium Bases



UniBase™ - Straight Ti-Base

H Gingival Height	0.5mm	1.5mm	3.0mm	0.5mm	1.5mm	3.0mm
Ø	5.2mm	5.2mm	5.2mm	5.2mm	5.2mm	5.2mm
Internal Hex	Anti-Rotational / Engaged			Rotational / Non-Engaged		
Zimax® WP	TBC-SV-WP	TBC-SV-WP15	TBC-SV-WP30	TBC-SV-WPR	TBC-SV-WP15R	TBC-SV-WP30R

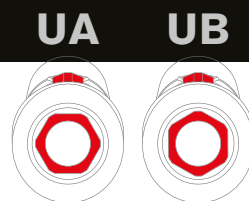


UniBase™ System

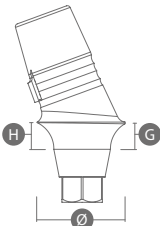
By Nikola Khoury

20° Angled Ti-Bases

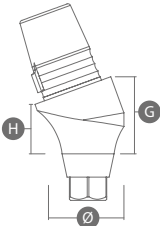
Patent Pending

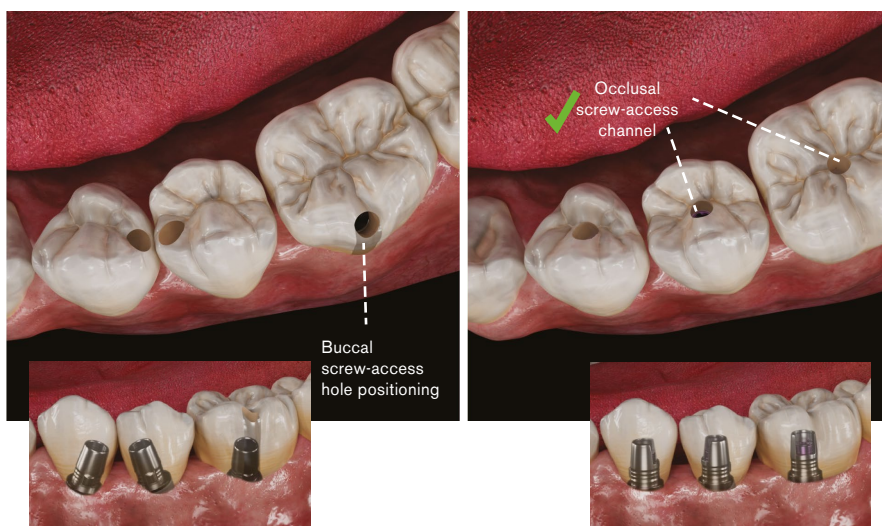
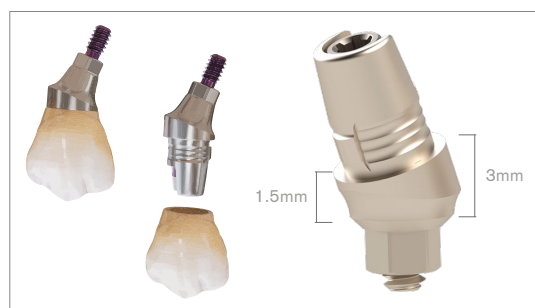


20° Angled Ti-Base - Anterior (U Type)

	UA	UB
		
H / G	0.7mm/0.7mm	0.7mm/0.7mm
Ø	4.5mm	4.5mm
Internal Hex		
Zimax® WP	TBC-SV-WP2007UA	TBC-SV-WP2007UB

20° Angled Ti-Base - Standard (A Type)

	AA	AB
		
H / G	1.5mm/3.0mm	1.5mm/3.0mm
Ø	5.2mm	5.2mm
Internal Hex		
Zimax® WP	TBC-SV-WP2015A	TBC-SV-WP2015B



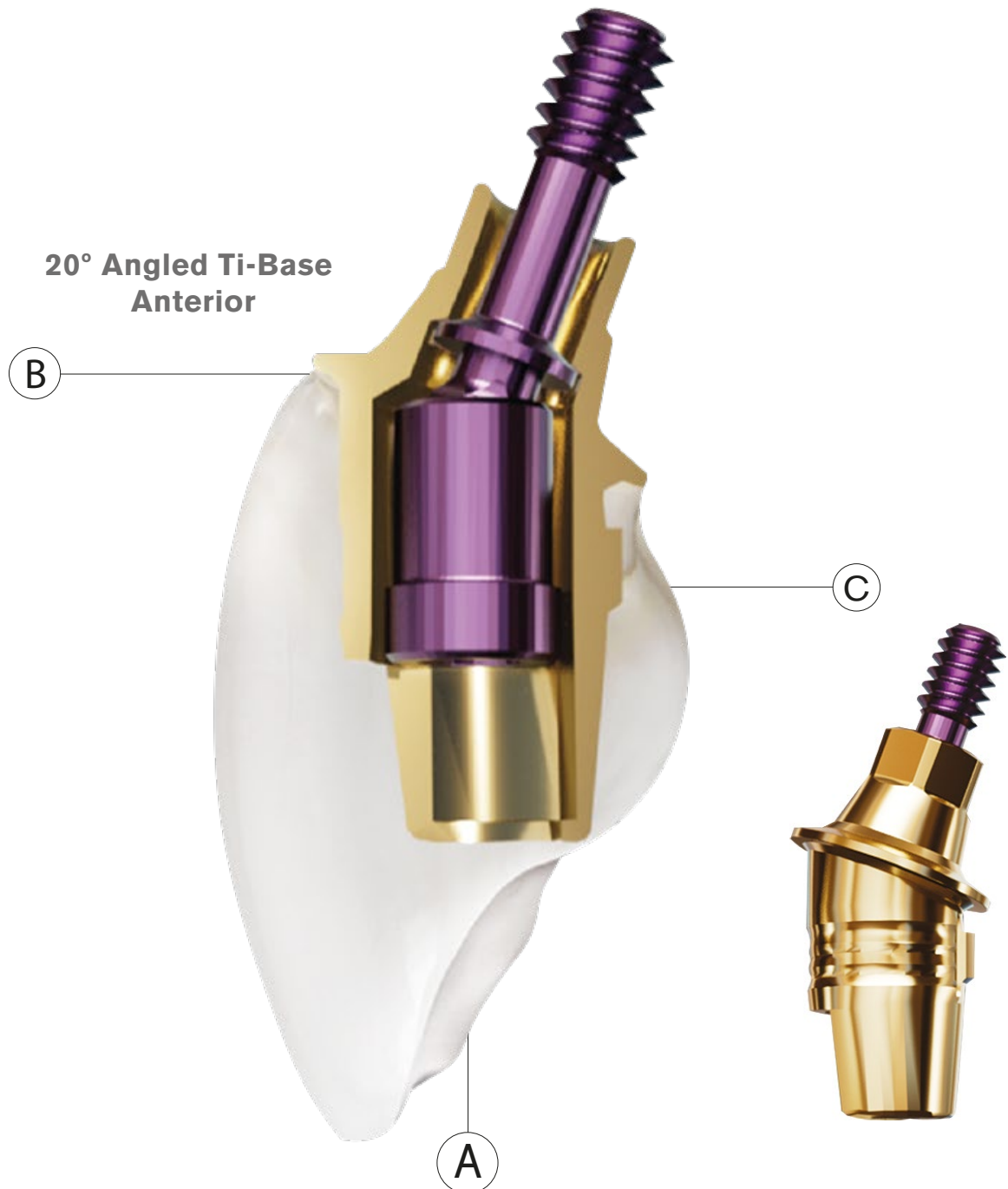
UniBase™ System

By Nikola Khoury

20° Angled Ti-Bases - Characteristics



⚠ Max. Tightening
Torque 22-25 NCM



**20° Angled Ti-Base
Anterior**

B

C

A

A-Simplicity!

A unique architectural design allows for the use of a standard screwdriver when mounting the UniBase, no special tools necessary.

B-Engineering!

The flat surface of the shoulder makes for an angulated base that interacts with the implant system in the same manner as a non-angulated base, allowing laboratories to use standard CAD/CAM procedures during the milling process.

C-Innovation!

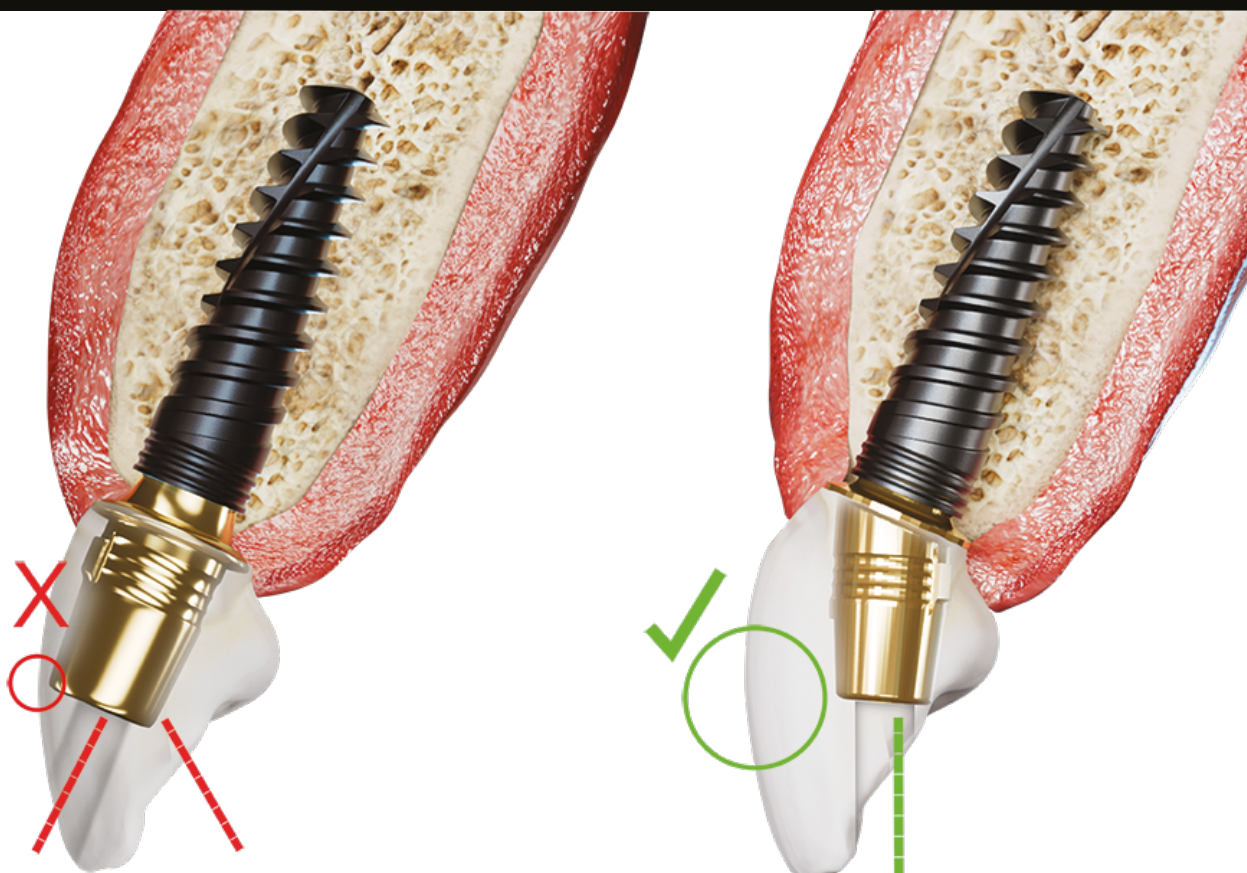
Our UniBase features an exclusive mechanism with a built-in "joint-shaft" that is secured into place during the fabrication process.

UniBase™ System

By Nikola Khoury

20° Angled Ti-Bases

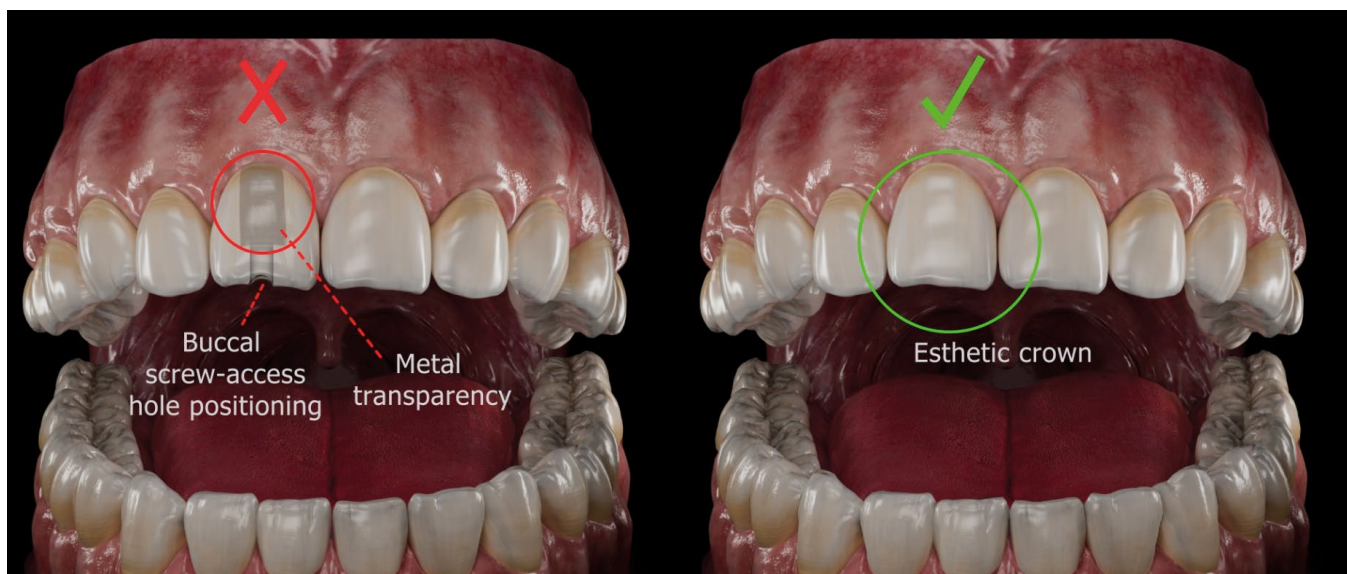
Patent Pending



Unfavorable buccal
screw-access channel

Lingual angled screw-ac-
cess channel

Ideal lingual
screw-access channel



Buccal
screw-access
hole positioning

Metal
transparency

Esthetic crown

Lab Analogs

Digital Implant Analogs



Digital Implant Analog

L	10mm
Ø	4.5mm
Internal Hex	
Zimax® WP	DIA-SV-WP

Connection Chart / Codes

Conical Connection	Connection code
ARTIC® NP	NA-NP
ARTIC® RP	NA-RP
FORTIS® RP	TS-RP
iOS3® SP	OS-SP
MORSY® RP	AR-RP
T-Oss® NP/Mini	TS-NP
Astix® WP	OS-WP
Nyo® RP	GM-RP
Gross® NC	CF-NC
Gross® RC	CF-RC
Internal Hex	Connection code
LAMINA® RP	SV-RP
Lgc® NL	SV-NL
Mass® NP	SV-NP
Zimax® WP	SV-WP



* Digital Analogs

NP / NC / NL / SP Platform Ø 3.5mm
 RP Platform Ø 4.0mm
 WP Platform Ø 4.5mm

* Straight TiBases

NP / NC / NL / SP Platform Ø 4.2mm / 4.5mm
 RP Platform Ø 4.2mm / 4.5mm
 WP Platform Ø 4.8mm / 5.2mm

* Angled TiBases

NP / NC / NL / SP Platform Ø 4.5mm
 RP Platform Ø 4.5mm
 WP Platform Ø 4.8mm / 5.2mm

Enjoy our
Multi Unit
System

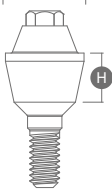


Multi Unit System

Classic Line M1.4



Straight MUA Abutment

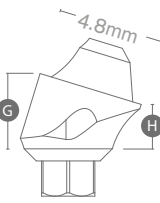
						
H	1mm	2mm	3mm	4mm	5mm	6mm
Internal Hex						
Zimax® WP	MUL-SV-WPNB10	MUL-SV-WPNB20	MUL-SV-WPNB30	MUL-SV-WPNB40	MUL-SV-WPNB50	MUL-SV-WPNB60

M1.4
Fixation Screw
Max. Torque 15NCM



17° Angled - Classic Line M1.4

30° Angled

						
H/G	1mm/2.5mm	2mm/3.5mm	3mm/4.5mm	1mm/3.5mm	2mm/4.5mm	3mm/5.5mm
Internal Hex	17°	17°	17°	30°	30°	30°
Zimax® WPZ	MUL-SV-WPNB1710	MUL-SV-WPNB1720	MUL-SV-WPNB1730	MUL-SV-WPNB3010	MUL-SV-WPNB3020	MUL-SV-WPNB3030

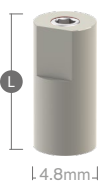
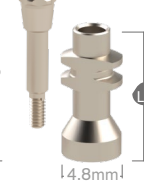
Healing Cap

MUA Scan Body

Closed Transfer

Open Transfer

Temporary Sleeve

				
L 5.2mm	L 10mm	L 7.2mm	L1/L2 10.7/4.2mm	L 12mm
Hex 1.25/Unigrip	Hex 1.25/Unigrip	Hex 1.25/Unigrip	Hex 1.25/Unigrip	Hex 1.25/Unigrip
MUL-NB-HAUG	MUL-NB-SBUG	MUL-NB-PUPUG	MUL-NB-OTUG	MUL-NB-SLVTUG

Digital Analog

Titanium Base

Adhesive Cap

Fixation Screw

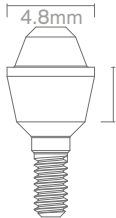
NikolasBase™

				
L 12mm	L 4.8mm	L 7.5mm	L 4.6mm	L 5.2mm
REF#	Hex 1.25/Unigrip	Hex 1.25/Unigrip	Hex 1.25/Unigrip	Hex 1.25/Unigrip
MUL-NB-DAN	MUL-NB-TBRUG	MUL-NB-ADRUG	MUL-NB-TSUG	MUL-NB-TBUG

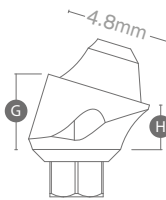
Multi Unit System

Multi-Clcq™ Line M1.6

Straight MUA Abutment - Multi-Clcq™ Line M1.6

						
H	1mm	2mm	3mm	4mm	5mm	6mm
Internal Hex						
Zimax® WP	MUL-SV-WPKS10	MUL-SV-WPKS20	MUL-SV-WPKS30	MUL-SV-WPKS40	MUL-SV-WPKS50	MUL-SV-WPKS60

17° Angled - Multi-Clcq™ Line M1.6

		
H/G	1mm/2.5mm	2mm/3.5mm
Internal Hex	17°	17°
Zimax® WP	MUL-SV-WPKS1710	MUL-SV-WPKS1720

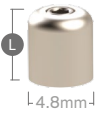
30° Angled

	
1mm/3.5mm	2mm/4.5mm
30°	30°
MUL-SV-WPKS3010	MUL-SV-WPKS3020

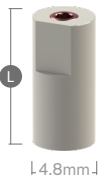
M1.6
Fixation Screw
Max. Torque 22NCM



Healing Cap


L 5.2mm
Hex 1.25/Unigrip
MUL-KS-HAUG

MUA Scan Body


L 10mm
Hex 1.25/Unigrip
MUL-KS-SBBUG

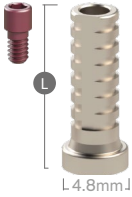
Closed Transfer


L 7.2mm
Hex 1.25/Unigrip
MUL-KS-PUPUG

Open Transfer


L1/L2 10.7/4.2mm
Hex 1.25/Unigrip
MUL-KS-OTUG

Temporary Sleeve


L 12mm
Hex 1.25/Unigrip
MUL-KS-SLVTIUG

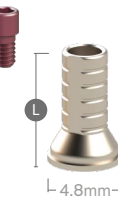
Digital Analog


L 12mm
REF#
MUL-KS-DAN

Titanium Base


L 4.8mm
Hex 1.25/Unigrip
MUL-KS-TBRUG

Adhesive Cap


L 7.5mm
Hex 1.25/Unigrip
MUL-KS-ADRUG

Fixation Screw

Max. 15NCM

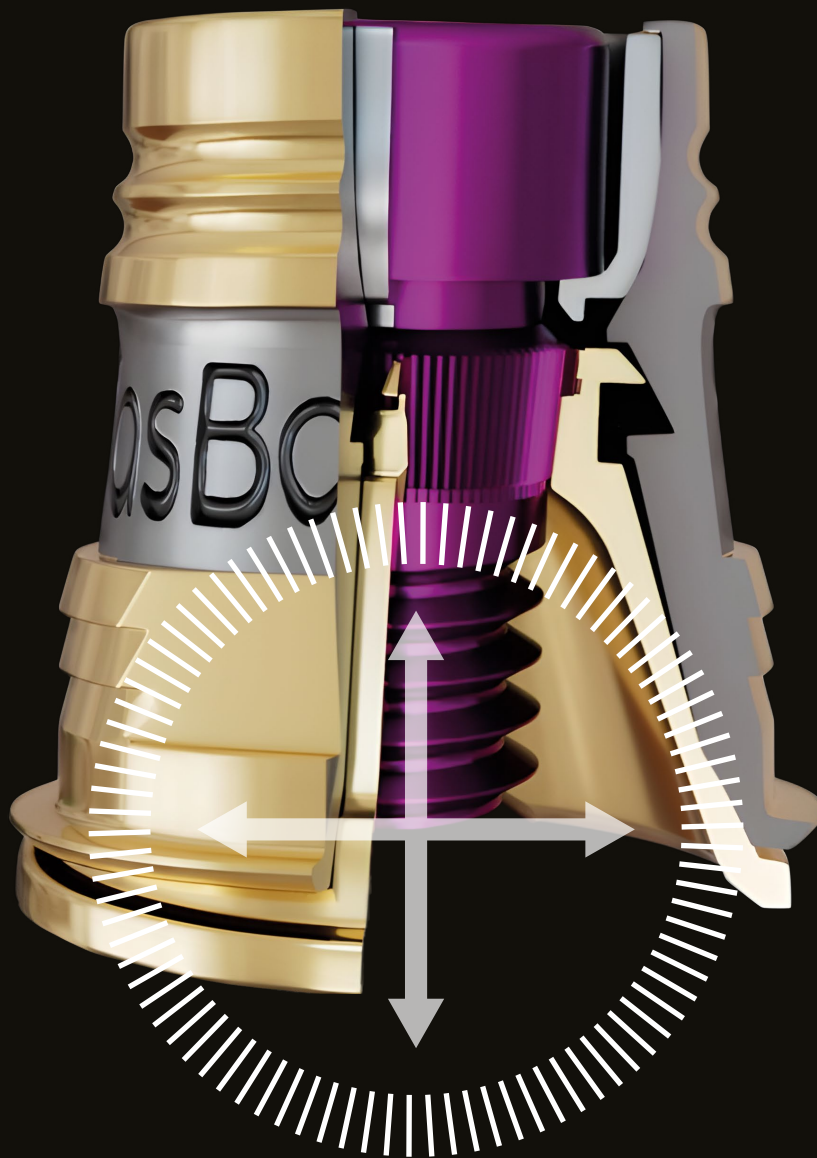
L 4.6mm
Hex 1.25/Unigrip
MUL-KS-TSUG

NikolasBase™

Max. 25NCM

L 5.2mm
Hex 1.25/Unigrip
MUL-KS-TBUG

NikolasBase™



Self-Adjustable Adhesive Cap

NikolasBase™

By Nikola Khoury

Self-Adjustable MUA Adhesive Cap



Standard MUA Ti-Base / Adhesive Cap

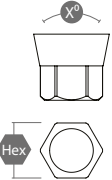
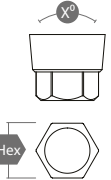
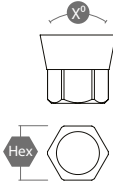
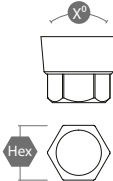
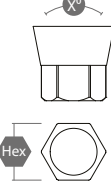
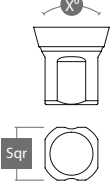
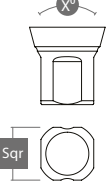


NikolasBase™ MUA Ti-Base / Adhesive Cap

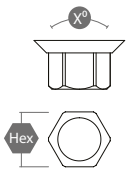
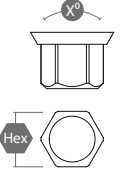
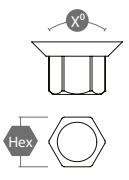
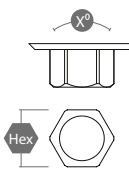


Platform Dimensions

Conical Connection

 X° - 24° Hex - 2.25mm ARTIC® NP NA-NP	 X° - 24° Hex - 2.65mm ARTIC® RP NA-RP	 X° - 22° Hex - 2.5mm FORTIS® RP TS-RP	 X° - 22° Hex - 2.1mm iOS3® SP OS-SP	 X° - 10° Hex - 2.35mm MORSY® RP AR-RP
 X° - 22° Hex - 2.1mm T-Oss® NP TS-NP	 X° - 16° Hex - 2.0mm Nyo® RP GM-RP	 X° - 22° Hex - 2.5mm Astix® WP OS-WP	 X° - 30° Sqr - 2.0mm Gnoss® NC CF-NC	 X° - 30° Sqr - 2.5mm Gnoss® RC CF-RC

Internal Hex

 X° - 90° Hex - 2.43mm LAMINA® RP SV-RP	 X° - 60° Hex - 2.0mm Lgc® NL SV-NL	 X° - 90° Hex - 2.1mm Mass® NP SV-NP	 X° - 90° Hex - 2.43mm Zimax® WP SV-WP
---	---	--	--

Screw Types

Recommended Tightening Torque

Screw Type

Hex 1.25mm



iOS3® SP
T-Oss® NP/Mini
Lgc® NL
Mass® NP
FORTIS® RP
Astix® WP
LAMINA® RP
Zimax® WP

Hex 1.25mm



Hex/Unigrip (Dual)

Unigrip **UG**



ARTIC® NP
ARTIC® RP

AR
 Hex 1.2mm



MORSY® RP

Torx



Gnoss® NC
Gnoss® RC

GM



Nyo® RP

Tightening Torque

Scan Bodies / Digital Transfers:

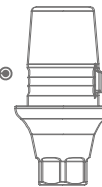
A/B/C 8-10 NCM



Straight Unibase™ Abutments:

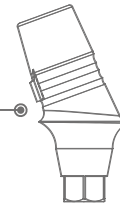
A 22 NCM

B/C 30 NCM



Angled Unibase™ Abutments:

A/B/C 22 NCM



Multi-Clicq™ - Multi-Unit Abutments:

A 22 NCM Straight/Angled

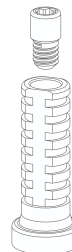
B/C 30 NCM Straight/Angled



Multi Unit Components:

Multi-Clicq M1.6 22 NCM

Classic Multi Unit M1.4 15 NCM



A

ARTIC® NP

iOS3® SP

T-Oss® NP

Nyo® RP

Gnoss® NC

Gnoss® RC

Lgc® NL

Mass® NP

B

ARTIC® RP

FORTIS® RP

MORSY® RP

Astix® WP

C

LAMINA® RP

Zimax® WP

a solution you
can **count** on.



 **Ritter Implants GmbH & Co. KG**

Freiburger Str. 45, 88400 Biberach, Germany
Tel: +49.7351.5292510

Headquarters: Edison Medical
4310 West Avenue, San Antonio, TX 78213, USA
Tel: +1.210.9793389 | edisonmed.com | zinon.dental